

power one ACCU *plus*
Mobile energy for optimum
hearing.

powerone)))



www.powerone-batteries.com

power one ACCU *plus* rechargeable hearing aid batteries.

Your demands

You want hearing aid batteries that you can rely on because of their outstanding performance:

power one ACCU *plus* are highly efficient rechargeable batteries based on NiMH technology characterized by the highest quality and reliability.

Your environmental contribution

The quickly-charging rechargeable batteries are also an ecological first choice. The use of rechargeable hearing aid batteries protects the environment.

power one rechargeable batteries are available in the sizes p13 accu, p312 accu, p10 accu and p675 accu and can be used in all hearing aids with the corresponding size.




power one ACCU *plus* fulfills your expectations:

- The stainless steel housing ensures high stability and safe corrosion protection.
- A quality feature is surface coating through nanotechnology – providing a unique, additional corrosion protection even with strong perspiration.
- Unlimited functionality under the most varied environmental conditions.
- Environmentally-friendly technology.
- 0% mercury, 0% lead, 0% cadmium.
- Fast-charging capabilities when used with the **power one pocketcharger**.
- No memory effect.
- Made in Germany.

Convenient packaging

power one ACCU plus is available in easy to handle packaging that stores your cells safely:

- The practical dispenser wheel facilitates the optimum removal of the rechargeable batteries.
- Maximum product protection thanks to the quality seal.
- The reclosable strap ensures safe storage.
- Good visible colour coding.
- Environmentally friendly PET and cardboard materials.

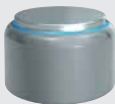

power one fast chargers for your mobility.

pocketcharger – ultra fast charging in credit card format

The **power one pocketcharger** system is suitable for charging one or two rechargeable **power one ACCU plus** hearing aid batteries.

Depending on which corresponding battery compartment is selected, rechargeable batteries of the sizes p13 accu, p312 accu or p10 accu can be charged. The **pocketcharger** is operated by a rechargeable lithium-polymer battery and therefore does not need to be plugged in. Universal power supply for recharging is included.

Depending on type, discharged cells are fully recharged within only 2.5 hours. A safety function, the refreshing protection, drains the still partially charged rechargeable batteries and thus improves the shelf life of the cells.



power one 675charger – a sturdy charging system

The **power one 675charger** serves to charge one or two **power one ACCU plus** type p675 accu rechargeable hearing aid batteries.

The discharged cells are completely recharged within 5 hours. For longer shelf life of the cells, the rechargeable batteries that are still partially charged are automatically discharged before the charging process.


The user-friendly **pocketcharger** and **675charger** fast chargers are **power one** products and therefore especially reliable and safe.



Please take note that:

The **power one** ACCU *plus* rechargeable batteries are recommended for the **pocketcharger** and **675charger**. Do not store ACCU *plus* in a discharged state.

Contact takes place through
top and bottom

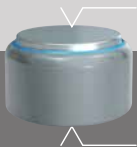

All **power one** rechargeable batteries are made in Ellwangen, Germany, in the world's largest and most advanced production facility and are subject to the toughest safety and quality control tests.

www.powerone-batteries.com