

# Phonak Naída™ Link

## Product information



Phonak Naída Link is a water resistant power hearing aid that works with the Advanced Bionics Naída CI Q90 or Q70. Since both the Phonak Naída Link and the Naída CI have Binaural VoiceStream Technology™, the client hears speech better in noise and on the telephone. Phonak Naída Link comes in a RIC and a UP model and offers the same choices of colors as Naída CI.

For more information visit [www.phonakpro.com](http://www.phonakpro.com)

# Product description



Phonak Naída Link RIC




Phonak Naída Link UP

		xS	xP	xSP+	HE7	HE7 680	Power SlimTube
Maximum Power Output (dB SPL)	2cc	109	123	130	139	132	139
	Ear simulator	121	132	137	144	139	144
Maximum gain (dB)	2cc	45	55	65	82	75	82
	Ear simulator	56	65	72	85	80	85
Frequency range (Hz)		<100 - 8900	<100 - 6200	<100 - 5100	<100 - 4900	<100 - 4900	<100 - 4900
Battery size		13	13	13	675	675	675
Push button (program / volume)		•	•	•	•	•	•
Volume control					•	•	•
Telecoil		•	•	•	•	•	•
IP rating: IP67 <sup>1</sup>		•	•	•	•	•	•
Direct Audio Input (DAI)		•	•	•	•	•	•
Hook					•	•	•
Power SlimTube					•	•	•
Secured battery compartment		•	•	•	•	•	•


# Fitting range

Phonak Naída Link RIC  
Standard xReceiver (xS)

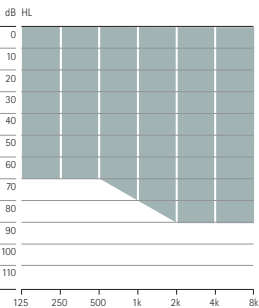


Mild to moderately-severe hearing loss, all audiometric configurations.

Phonak Naída Link RIC  
Power xReceiver (xP)

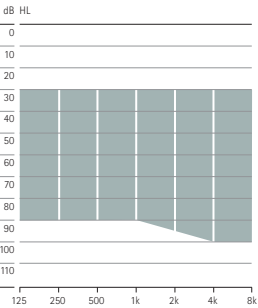


Mild to severe hearing loss, all audiometric configurations.

Phonak Naída Link RIC  
SuperPower plus xReceiver (xSP+)

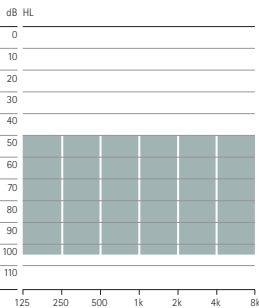

Moderate to profound hearing loss, all audiometric configurations.

Phonak Naída Link UP

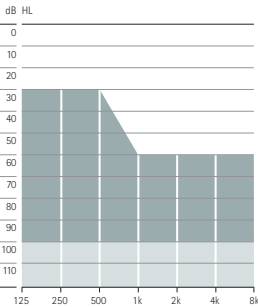

Severe to profound hearing loss, all audiometric configurations.

<sup>1</sup> IP67 indicates that the hearing aid is water resistant and dust tight. It survived immersion in 1m of water for 30 minutes and 8 hours in a dust chamber as per the IEC60529 standard. No traces of dust were evident within the housing.

# Housing colors



# Acoustic coupling










Hook / SlimTube	Hook	Available as damped (HE7 680), undamped (HE7) and as pediatric hook (mini hook HE7 680)	
	Power SlimTube	Available in left, right and five different lengths (00, 0, 1, 2, 3), requires custom earmold	
xReceiver	Length	Available in four different lengths (0, 1, 2, 3), left and right	
	Dome	Open dome	Available in three sizes (S, M, L)
		Closed dome	Available in three sizes (S, M, L)
		Power dome	Available in three sizes (S, M, L)
	SlimTip*	Hard material	Hollow shell
			Solid shell
		Soft material	Solid shell
	cShell*	Hard material	

\*Available with extra canal length option

# Fitting




Software	Phonak Target™ 4.3 or higher <ul style="list-style-type: none"><li>Consult the Bimodal Fitting Report from the Naida CI Fitting Software (SoundWave) for the fitting of Naida Link.</li><li>Follow the fitting instructions on the Phonak Naida Desktop Fitting Guide.</li></ul>
Interfaces	iCube II, NOAHLINK, HI-PRO, HI-PRO2

# Roger portfolio

Roger receivers	Phonak Naída Link RIC	Phonak Naída Link UP
 Roger 10		•
 Roger 15	•	
 Roger X / AS10		•
 Roger X / AS15	•	
 Roger MyLink	•	•
Roger microphones	Phonak Naída Link RIC	Phonak Naída Link UP
 Roger Pen	•	•
 Roger EasyPen	•	•
 Roger Clip-on Mic	•	•
 Roger Table Mic	•	•

For information on Roger microphones please visit [www.phonakpro.com/roger](http://www.phonakpro.com/roger)

# Wireless accessories

		Phonak Naída Link RIC	Phonak Naída Link UP
	Phonak DECT Phone	•	•
	Phonak ComPilot <sup>1</sup>	•	•
	Phonak TVLink II	•	•

<sup>1</sup> Phonak ComPilot needs to be fitted to the Naída CI Q90 or Q70 by using the SoundWave fitting software from Advanced Bionics.

## Feature overview

	Phonak Naída Link RIC	Phonak Naída Link UP
StereoZoom •	•	•
ZoomControl •	•	•
DuoPhone •	•	•
UltraZoom	Premium	Premium
SNR-Boost	•	•
SoundRecover	•	•
SoundFlow	5	5
Real Ear Sound	•	•
Channels	20	20
WhistleBlock	•	•
NoiseBlock	•	•
WindBlock	•	•
EchoBlock	•	•
SoundRelax	•	•
QuickSync (wireless)	•	•
Manual programs	5	5
Streaming programs (wireless)	5	5
User Preference Tuning	•	•
auto Acclimatization	•	•
Wireless fitting	•	•

• Binaural VoiceStream Technology™

---

# Life is on

We are sensitive to the needs of everyone who depends on our knowledge, ideas and care. And by creatively challenging the limits of technology, we develop innovations that help people hear, understand and experience more of life's rich soundscapes.

Interact freely. Communicate with confidence.  
Live without limit. Life is on.

[www.phonakpro.com](http://www.phonakpro.com)